

Knight®

by Midmark

Midmark Corporation • P O Box 286 • 60 Vista Drive
Versailles • Ohio • 45380-0286 • www.midmark.com
Phone: 937-526-3662 Fax: 937-526-5542

ESCORT FIXED SIDE MOUNT FOR KNIGHT CHAIRS

KIT NO. 151665 INSTALLATION INSTRUCTIONS

WARNING: ALTERING THE CENTER OF GRAVITY OF AN UNSECURED CHAIR MAY CAUSE THE CHAIR TO TIP. SECURE THE CHAIR TO THE FLOOR BEFORE THE DENTAL UNIT, OR ANY OTHER ACCESSORY TO THE CHAIR. PERSONAL INJURY MAY RESULT IF THE CHAIR IS NOT PROPERLY SECURED.

NOTE: WHEN INSTALLING THE UNIT FOR LEFT-HAND OPERATION, IT IS NECESSARY TO REPOSITION THE SUPPORT ASSY ON THE CHAIR MOUNT . SEE FIGURE 2.

1. RAISE THE CHAIRTOP TO A CONVENIENT WORKING HEIGHT. SECURE THE UNIT MOUNT TO THE CHAIR USING (2) HEX SOCKET HEAD SCREWS. TIGHTEN THE MOUNTING SCREWS SECURELY.

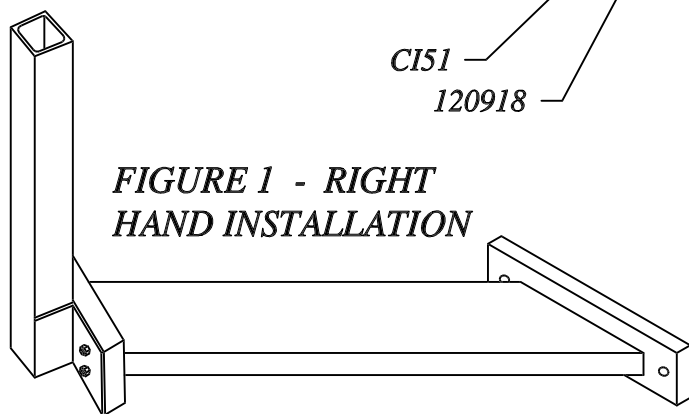


**FIGURE 1 - RIGHT
HAND INSTALLATION**

CI51

120918

**FIGURE 2 - LEFT
HAND INSTALLATION**

117990

2. ASSEMBLE THE EQUIPMENT TO THE MOUNT AS DIRECTED IN THE EQUIPMENTS MANUAL.

THE FIXED SIDE MOUNT KIT INCLUDES A UNIT MOUNT & HARDWARE

TP1164 Rev. B (10/99)